



SAFETY STORAGE CABINETS

CONTENTS

Safety Storage Cabinet Features 3

Model

Flammable Safety Cabinets	5
Corrosive Safety Cabinets	6
Combustible Safety Cabinets	7
Hazardous Safety Cabinets	8
Polypropylene Safety Cabinets	9
Stainless Steel Safety Cabinets	10
Flammable Drum Storage Cabinets	11


Warranty 12

Safety Storage Cabinet Features



- Constructed from industrial grade electro-galvanised steel 18 gauge for strength and durability
- Epoxy powder coated to meet harsh industrial/laboratory environments
- Double layered construction with 1.5" airspace between layers
- Fully welded 2" sill to contain chemical spills within the cabinet
- 2 ventilation ports protected with flash arrestors
- 3 point door latching mechanism
- Large Warning Labels for easy identification
- Adjustable shelves with 1/2" centers
- Complies and 3rd party lab tested to NFPA 30 6.3.3 para 2

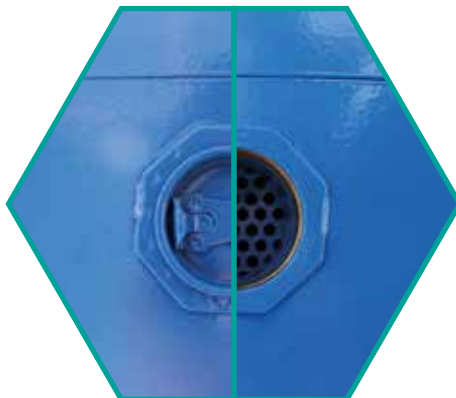
Adjustable galvanised steel shelves with heavy duty stainless steel pins



Double wall construction



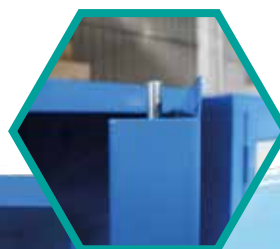
2" vents with flash arrestors (exterior)



3 pin lock



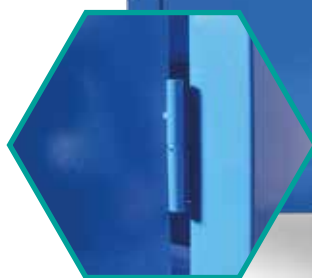2" deep spill sump



Heavy duty hinges

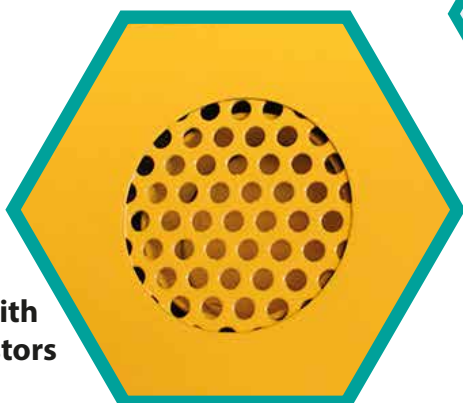


Optional PP tray






3 point door latch




2" vents with flash arrestors (interior)








Flammable Safety Cabinets

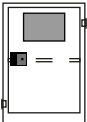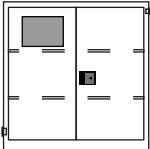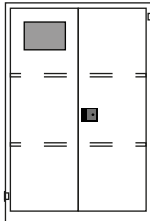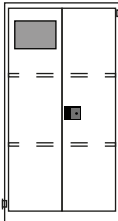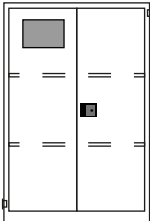
MODEL	AD15F	AD30F	AD45F	AD60F	AD90F
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

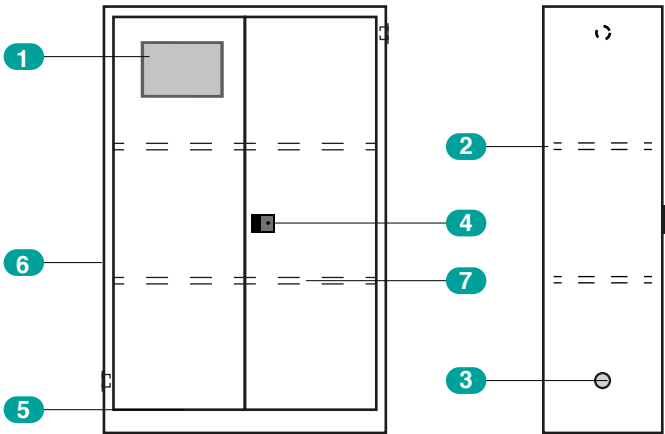

1. Safety Warning Sticker
2. Double wall construction with 1-1/2" air space
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. Epoxy powder coated
7. Galvanised steel shelves adjusted on 1/2" centers

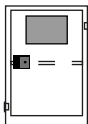
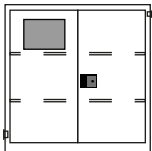
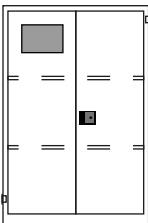
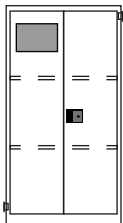
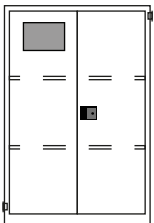
Corrosive Safety Cabinets


MODEL	AD15C	AD30C	AD45C	AD60C	AD90C
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

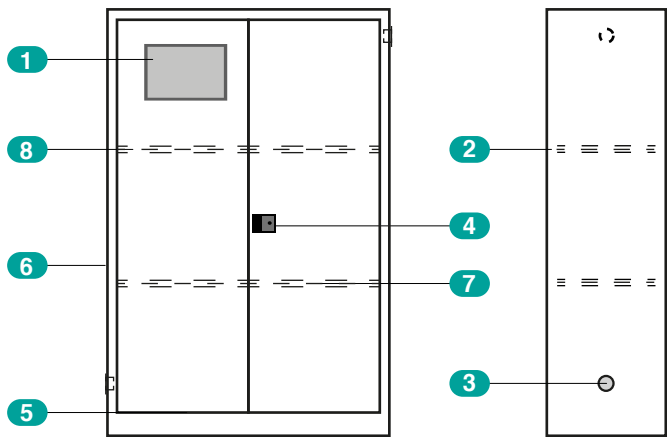

1. Safety Warning Sticker
2. Double wall construction with 1-1/2" air space
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. Epoxy powder coated
7. Galvanised steel shelves adjusted on 1/2" centers
8. PP Trays

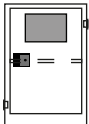
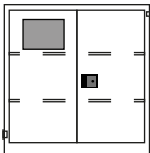
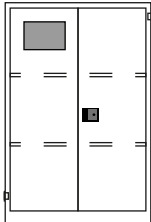
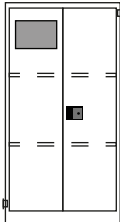
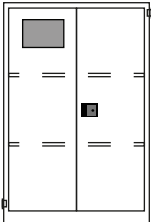
Combustible Safety Cabinets



MODEL	AD15R	AD30R	AD45R	AD60R	AD90R
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

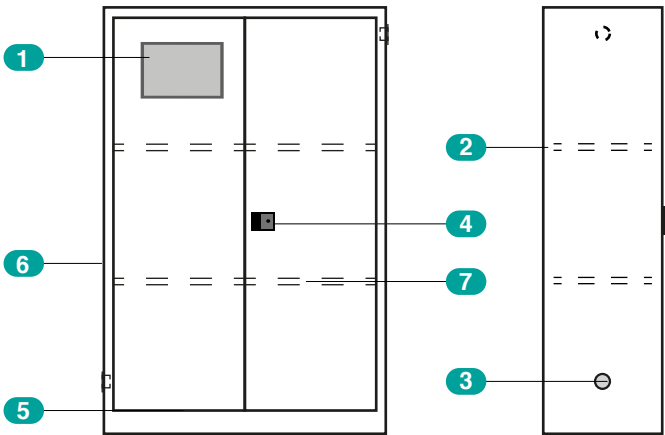

- 1. Safety Warning Sticker
- 2. Double wall construction with 1-1/2" air space
- 3. 2 nos 2" vents with flash arrestors
- 4. Flush mounted locking handle
- 5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
- 6. Epoxy powder coated
- 7. Galvanised steel shelves adjusted on 1/2" centers

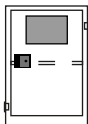
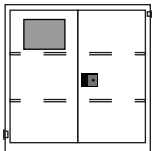
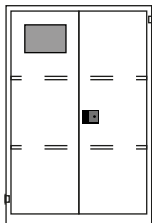
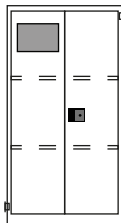
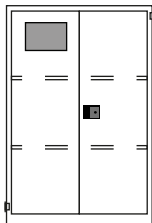
Hazardous Safety Cabinets


MODEL	AD15H	AD30H	AD45H	AD60H	AD90H
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

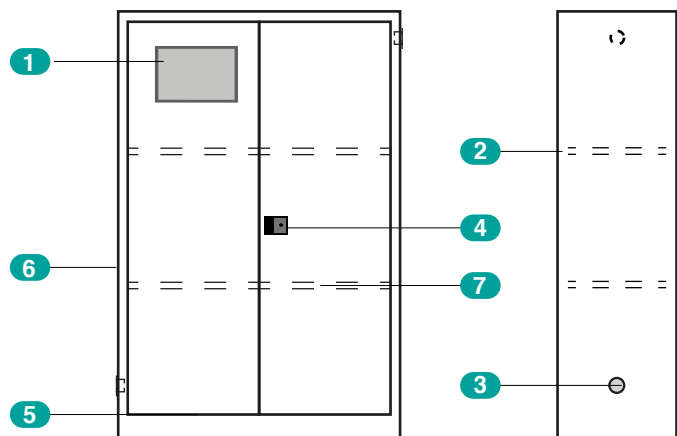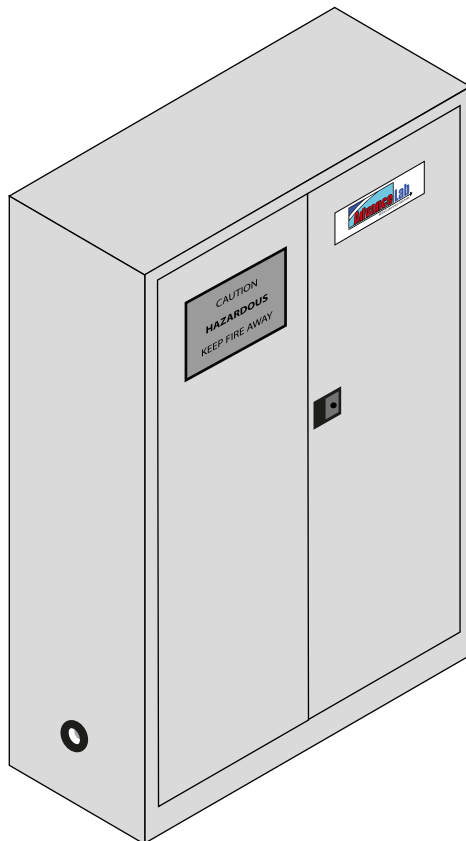

1. Safety Warning Sticker
2. Double wall construction with 1-1/2" air space
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. Epoxy powder coated
7. Galvanised steel shelves adjusted on 1/2" centers

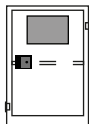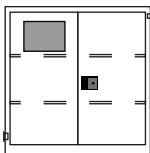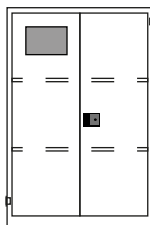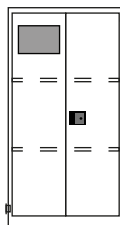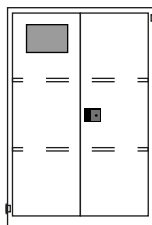Polypropylene Safety Cabinets






MODEL	AD15PP	AD30PP	AD45PP	AD60PP	AD90PP
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

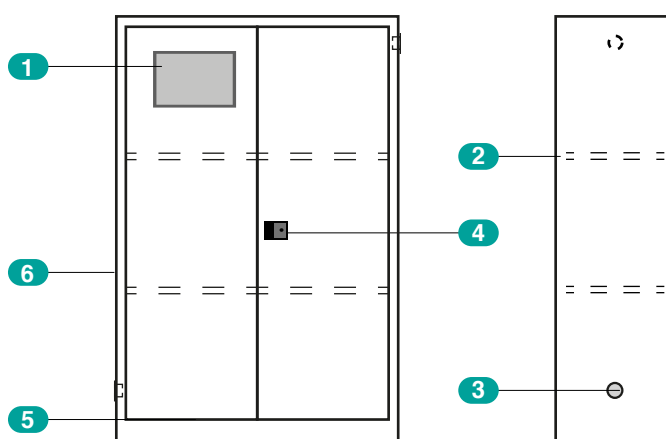


1. Safety Warning Sticker
2. Full PP hot weld construction (Single Layer)
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. PP adjustable shelves

* Single layer and not compliant to NFPA 30

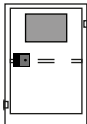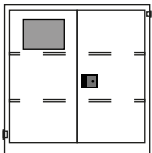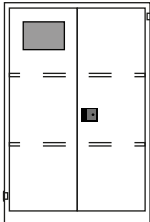
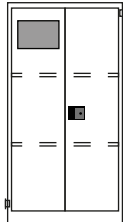
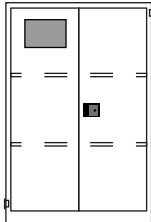
Stainless Steel Safety Cabinets



MODEL	AD15SS	AD30SS	AD45SS	AD60SS	AD90SS
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

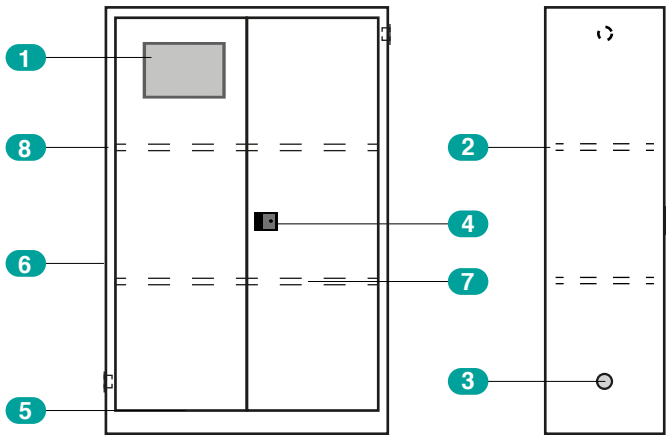

1. Safety Warning Sticker
2. Double wall construction with 1-1/2" air space
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. Stainless steel construction
7. Stainless steel shelves adjusted on 1/2" centers
8. Constructed from SS304 with 4B finish (SS316 Optional)

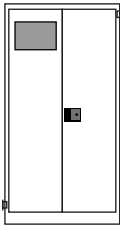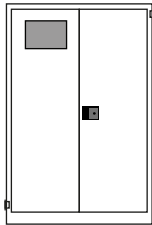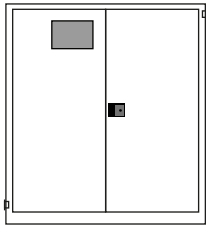Flammable Drum Storage Cabinets

MODEL	AD55DS1	AD30DS2	AD55DS2
Dimensions, mm (LxWxH)	864 x 864 x 1651 	1092 x 864 x 1651 	1496 x 864 x 1651 
Weight	150 Kg	220 Kg	280 Kg
Design	Manual closing door		
Door	2 door	2 door	2 door
Capacity	1 x 55 gallon drum	2 x 30 gallon drum	2 x 55 gallon drum

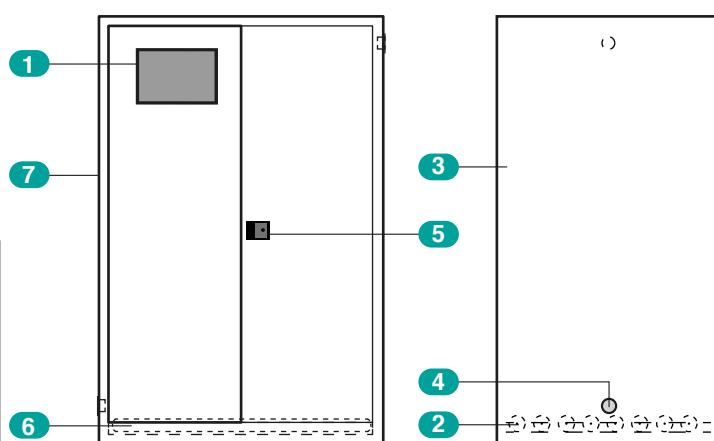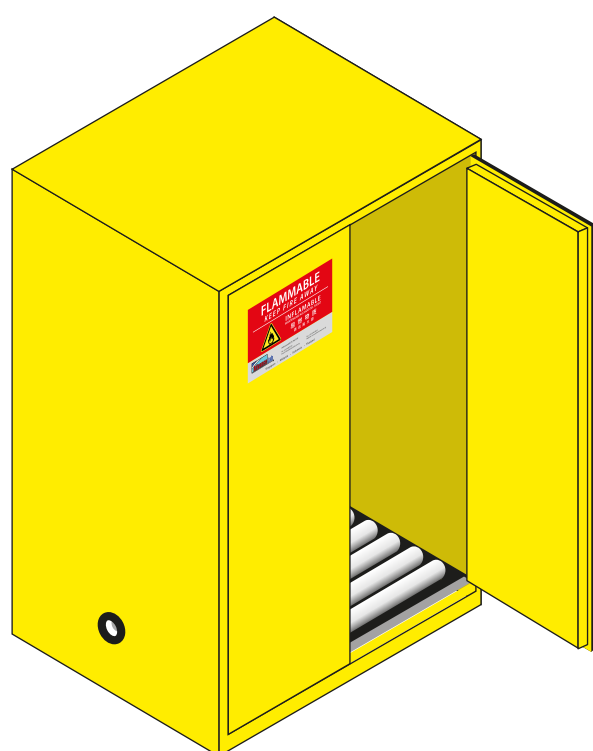

1. Safety Warning Sticker
2. 60mm dia rollers, steel
3. Double wall construction with 1-1/2" air space
4. 2 nos 2" vents with flash arrestors
5. Flush mounted locking handle
6. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
7. Epoxy powder coated

Warranty

Products Manufactured by Advancelab: Advancelab(S) Pte.Ltd., warrants products that it manufactures to be free from defects for a period of 12 months for parts, commencing from the date of shipment. Advancelab's sole responsibility is to repair or replace, at its option, any part of the product that proves defective or malfunctioning during this time limit. This warranty is void if the equipment is abused or modified by the customer, is operated outside Advancelab's operating instructions or specifications, or is used in any application other than that for which it is specified. This warranty does not include routine maintenance or service procedures, breakage, shipping damage, nor damage from misuse, intentional or unintentional abuse, neglect, natural disasters, or acts of God.

Freight Shortage or Damage: Upon receipt of any equipment from Advancelab, customer shall immediately unpack and inspect for damage or shortage. The customer shall not accept a damaged package or a short shipment until the carrier makes a "damage or shortage" notation on both the carrier's and customer's copy of the freight bill or delivery receipt. Service title passes when the shipment is loaded, so customer is responsible for filing and collecting a freight claim. Any replacement products must be ordered and paid for separately.

Generally, customers can improve the chance of collecting on a freight claim by following these procedures:

1. Formally requesting that the carrier inspect the shipment immediately upon suspecting damage or shortage to verify condition.
2. Notifying the carrier upon discovery of concealed damage and requesting an inspection within 15 days of receipt, both in person or phone and following up via mail.
3. Keeping the shipment as intact as possible, including retaining original packaging materials and keeping the product as close to the original receiving location as possible.
4. Holding salvage for disposition by the carrier.

All Claims: Advancelab (S) Pte.Ltd., expressly disclaims all other warranties, expressed or implied or implied by statute, including the warranties of merchantability or fitness for intended use. Advancelab is not responsible for consequential or incidental damages arising out of the purchase or use of the products supplied by Advancelab. Advancelab is not liable for damage to facilities, other equipment, products, property or personnel of others, or of their agents, suppliers, or affiliated parties, which is caused or alleged to have been caused by products supplied by Advancelab. In any event or series of events, Advancelab's total liability for any and all damages whatsoever is limited to the lesser of the actual damages or the original invoice cost of the items alleged to have caused the damage. The customer's sole and exclusive remedy for any cause of action whatsoever is repair or replacement of the non-conforming products or refund of the actual purchase price, at the sole option of Advancelab. All claims must be made in writing within 90 days of the date the product was shipped. Any claims not made within this time limit shall be deemed waived by the customer. Advancelab is not responsible for any additional costs of repair caused by poor packaging or in-shipment damage during return.

Warranty Returns: All warranty returns must be authorized in advance by Advancelab and approved by writing. Unless approved in advance for good reason, all returns must be in original condition, including all manuals, and must be packaged in original packaging materials. All returned goods are to be shipped to Advancelab, freight prepaid at customer's expense.



SAFETY STORAGE CABINETS

Singapore

Advancelab (S) Pte Ltd

253 Kaki Bukit Ave 1
Singapore 416061
Tel: +65 6448 8255
Fax: +65 6448 9833
Email: info@advancelab.com.sg

Malaysia

Advancelab Sdn Bhd

Johor Factory: No. 3388, Jalan
Pekeliling Tanjung 27/2, Kawasan
Perindustrian Indahpura, 81000
Kulaijaya, Johor, Malaysia.
Tel: +607 660 8877
Fax: +607 660 8866
Email: info-my@advancelab-global.com

Myanmar

Advancelab Scientific &
Engineering Co., Ltd

No.(81/1)-1A, Myin Thar 7th Street,
(14/1)Ward, South Okkalapa Township,
Yangon, Myanmar
Tel: +95 (9) 779753802
Email: info-mm@advancelab-global.com

Indonesia

PT. Advancelab Saintifik

253 Kaki Bukit Ave 1
Singapore 416061
Tel: +65 6448 8255
Fax: +65 6448 9833
Email: info@advancelab.com.sg

Thailand

Advancelab (Thailand)
Co., Ltd.

147 / 170 - 171 Baromrajchonnane
Road, Arunamarin Bangkoknoi, Bangkok
10700, Thailand.
Tel: +66 2 884 9480
Fax: +66 2 884 6441
Email: primasci@primasci.com

U.A.E (Dubai)

Advancelab FZCO

253 Kaki Bukit Ave 1
Singapore 416061
Tel: +65 6448 8255
Fax: +65 6448 9833
Email: info@advancelab.com.sg

Vietnam

M&T International
Trading Services
Co., Ltd

46A Phan Dinh Phung, Tan Thanh Ward,
Tan Phu District, Ho Chi Minh City
700000, Vietnam.
Tel: +84 28 6656 0610
Email: info@advancelab.com.sg

Website:

www.advancelabglobal.com



Singapore • Malaysia • Myanmar • Indonesia • Thailand • U.A.E • Vietnam