

WALK IN

# FUME HOOD



*biz*SAFE  
ST~~R~~

ACS  
REGISTRARS  
Health & Safety  
ISO 45001

ACS  
REGISTRARS  
Quality Assurance  
ISO 9001



<b>Standard Hood Features</b>	<b>3</b>
<b>Laboratory Fume Handling</b>	<b>4</b>
<b>Walk In Fume Hood -AD Series</b>	<b>5</b>
<b>Models</b>	<b>6</b>
<b>Technical Specification</b>	<b>7</b>
<b>Options and Accessories</b>	<b>8</b>
<b>Warranty</b>	<b>9</b>

# Standard Hood Features

- Double layered design
- Aerodynamic front fascia
- Industrial grade galvanised steel construction
- Heavy epoxy powder coated finish to withstand harsh laboratory environments
- Phenolic resin inner liner and baffle with built in slots for efficient scavenging
- Purpose built light box external to airflow, preventing any corrosion and minimizing turbulence
- Electronic light ballast for flicker free usage and user eye comfort
- Counterbalanced sash assembly for years of trouble free use
- Tempered safety glass with Easygrip® hand rail
- 2 side access panels built into inner liner for future addition of utilities
- All fume hoods comply to exhaust requirements and face velocity test as per ASHRAE 110-95; when ordered with control systems, units also comply to BS 7258 and EN 14175

# Laboratory Fume Handling



Laboratory fume hoods or fume cupboards represent the primary line of defense for personnel against exposure to harmful, toxic and dangerous fumes or chemicals as they are produced.

**“Fume hoods are essentially aerodynamic scavenging devices that remove dangerous fumes from a laboratory and expel it into the environment where it is sufficiently diluted to no longer pose a health risk.”**

Advancelab Laboratory fume hoods have been manufactured to provide the highest levels of safety for laboratory personnel.

Fume hoods when sized with an appropriate chemical resistant centrifugal exhaust fan will protect laboratory staff from exposure to harmful fumes.

Please contact our distributors or offices to help or design an efficient exhaust system.

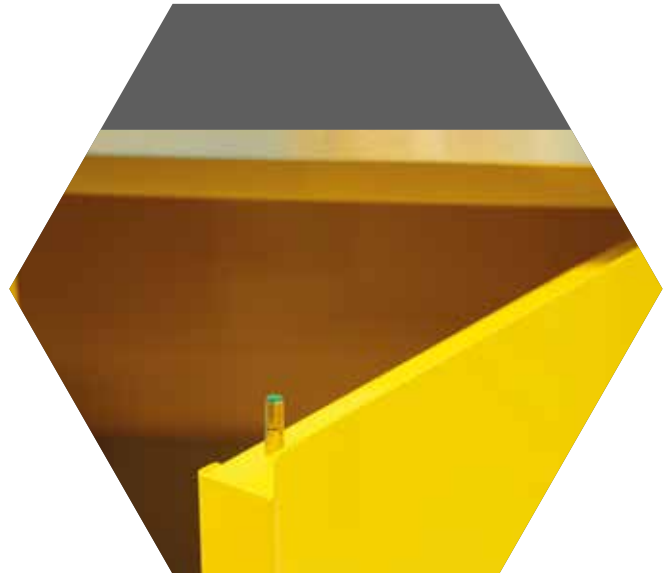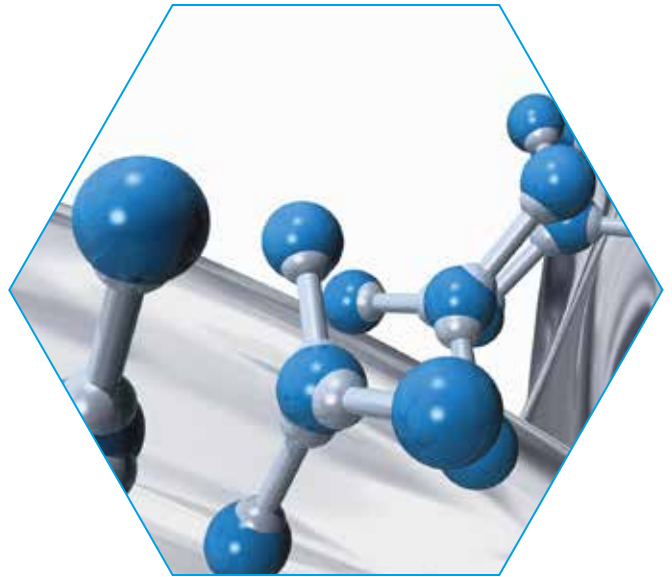

# Walk In Fume Hood

## AD Series


Walk In Fume Hoods




# Models

MODEL	AD04FWH	AD05FWH	AD06FWH	AD08FWH
Dimensions, (LxWxH)	1220 mm x 920mm x 2420mm	1520mm x 920mm x 2420mm	1820mm x 920mm x 2420mm	2420mm x 920mm x 2420mm
Body	Electrogalvanised steel with heavy epoxy powder coated finish			
Worktop	N/A			
Inner Liner & Baffle	Phenolic resin, laboratory grade			
Airfoil	N/A			
Standards Complied	ASHRAE 110-95 BS 7258 (with Optional Controller) EN 14175 (with Optional Controller)			
Services	<ul style="list-style-type: none"> <li>• 1 Pre Plumbed Cold Water Remote</li> <li>• 1 Pre plumbed Gas Remote</li> <li>• 4 x Pre wired 13 Amp Switch socket outlets (BS)</li> <li>• 1 PP Cup Sink</li> </ul>			
Controls	Fan ON/OFF switch Light ON/OFF switch			
Base	N/A			
Duct Collar	200 mm	250 mm	300 mm	300 mm x 2
Design Velocity	0.5 m/s at 500mm sash height			
Extract Volume at 500 mm Height Sash Opening	1188 CMH	1512 CMH	1836 CMH	2481 CMH


## Sash Configuration

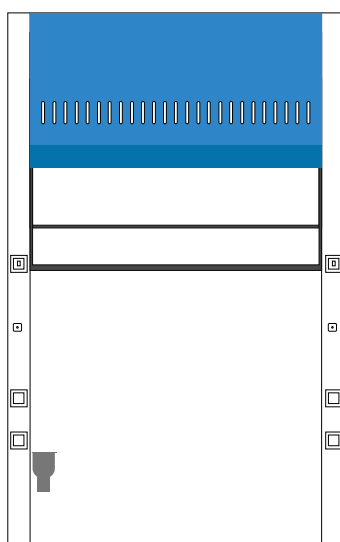

Working Position

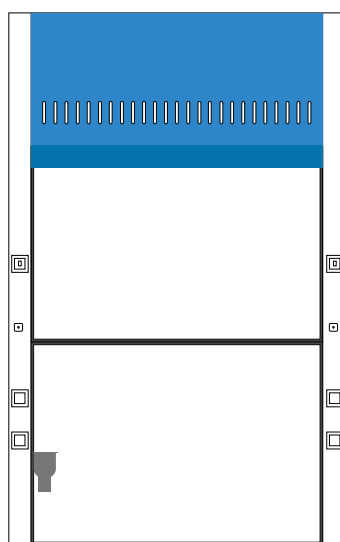


Full Close Position

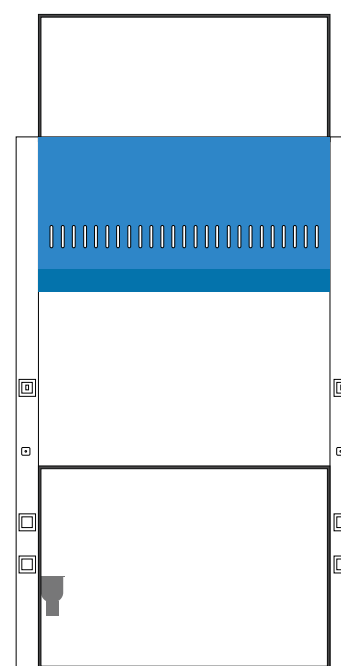


Full Open Position

# Technical Specifications

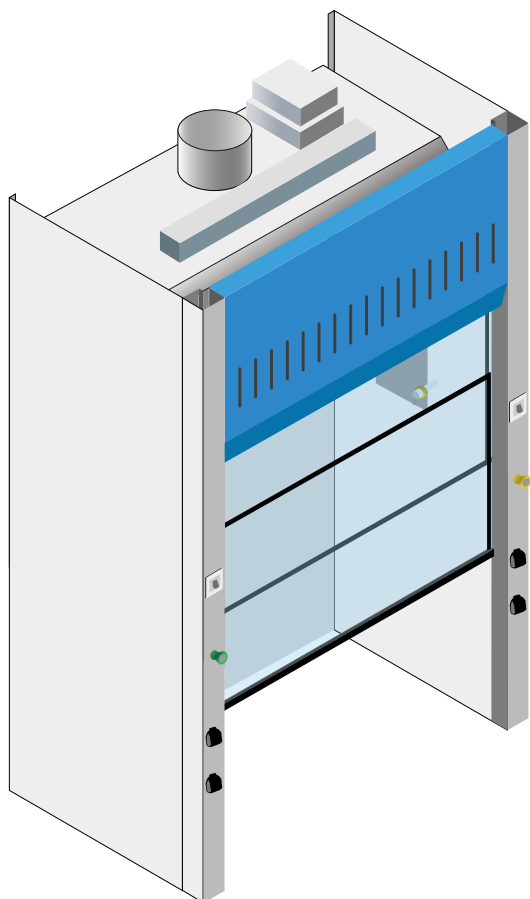


1. PVC Exhaust Collar
2. Electrical Junction Box
3. External Light Box
4. Side Access Panel
5. Removable Front Panel
6. Tempered Glass Sash
7. PP Cup Sink
8. Electrical Outlets
9. Fan ON/OFF Switch
10. Light ON/OFF Switch
11. Gas Fixture Remote
12. Water Fixture Remote
13. Gas Faucet
14. Water Faucet


Top View

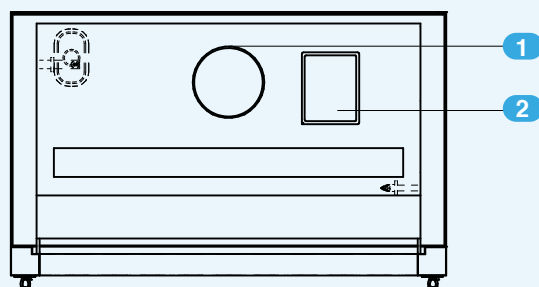

Side View

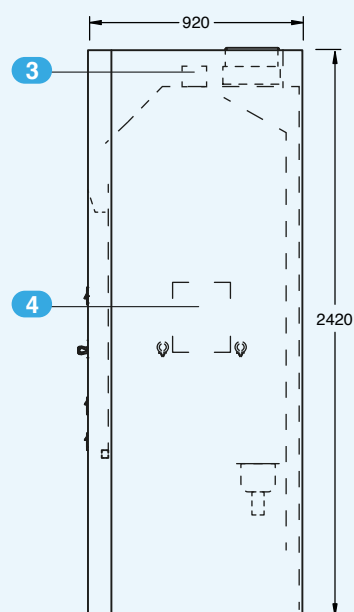

Front View

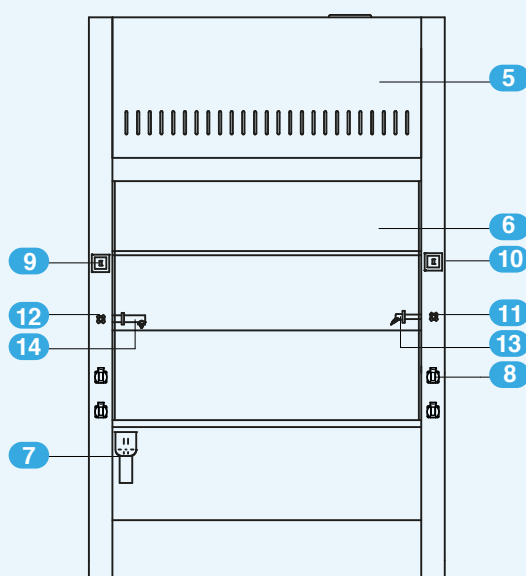

# Options and Accessories

## Superstructure

- Industrial Grade Electrogalvanised Steel (Standard)
- SS 304
- SS 316
- Polypropylene
- PVDF

## Distillation Grid

- Aluminium
- SS 304
- SS 316

## Inner Liner Construction


- Phenolic Resin (Standard)
- SS 304
- SS 316
- Polypropylene
- PVDF

## Sink

- Drip Cup
  - Polypropylene (Standard)
- Custom Sink
  - SS 304
  - SS 316
  - Polypropylene (Hot Weld Construction)
  - Size: \_\_\_\_\_ mm (L x W x H)  
*please specify*



- Location  
x = \_\_\_\_\_ mm  
y = \_\_\_\_\_ mm

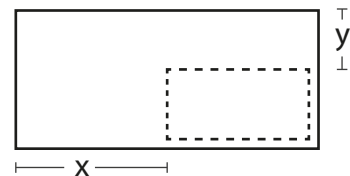

## Microprocessor Airflow Alarm

- Digital Display with Alarm
- Digital Display with Alarm & VAV/VSD Output



## Services

- Compressed Air (CA)
- Distilled Water (DI)
- Nitrogen (N2)
- Oxygen (O2)
- Vacuum (Vac)


## Ex Proof Configurations

- Ex proof lights
- Ex proof on/off switch

Please advice factory of details required

- Ex proof Electrical Sockets
- Ex proof electrical junction box

## Electrical Outlet

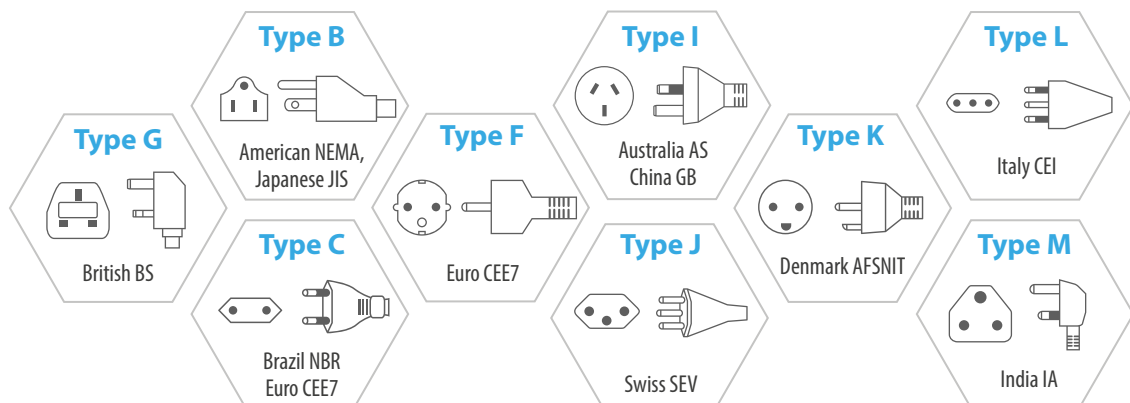



# Warranty

**Products Manufactured by Advancelab:** Advancelab(S) Pte.Ltd., warrants products that it manufactures to be free from defects for a period of 12 months for parts, commencing from the date of shipment. Advancelab's sole responsibility is to repair or replace, at its option, any part of the product that proves defective or malfunctioning during this time limit. This warranty is void if the equipment is abused or modified by the customer, is operated outside Advancelab's operating instructions or specifications, or is used in any application other than that for which it is specified. This warranty does not include routine maintenance or service procedures, breakage, shipping damage, nor damage from misuse, intentional or unintentional abuse, neglect, natural disasters, or acts of God.

**Freight Shortage or Damage:** Upon receipt of any equipment from Advancelab, customer shall immediately unpack and inspect for damage or shortage. The customer shall not accept a damaged package or a short shipment until the carrier makes a "damage or shortage" notation on both the carrier's and customer's copy of the freight bill or delivery receipt. Service title passes when the shipment is loaded, so customer is responsible for filing and collecting a freight claim. Any replacement products must be ordered and paid for separately.

Generally, customers can improve the chance of collecting on a freight claim by following these procedures:

1. Formally requesting that the carrier inspect the shipment immediately upon suspecting damage or shortage to verify condition.
2. Notifying the carrier upon discovery of concealed damage and requesting an inspection within 15 days of receipt, both in person or phone and following up via mail.
3. Keeping the shipment as intact as possible, including retaining original packaging materials and keeping the product as close to the original receiving location as possible.
4. Holding salvage for disposition by the carrier.

**All Claims:** Advancelab (S) Pte.Ltd., expressly disclaims all other warranties, expressed or implied or implied by statute, including the warranties of merchantability or fitness for intended use. Advancelab is not responsible for consequential or incidental damages arising out of the purchase or use of the products supplied by Advancelab. Advancelab is not liable for damage to facilities, other equipment, products, property or personnel of others, or of their agents, suppliers, or affiliated parties, which is caused or alleged to have been caused by products supplied by Advancelab. In any event or series of events, Advancelab's total liability for any and all damages whatsoever is limited to the lesser of the actual damages or the original invoice cost of the items alleged to have caused the damage. The customer's sole and exclusive remedy for any cause of action whatsoever is repair or replacement of the non-conforming products or refund of the actual purchase price, at the sole option of Advancelab. All claims must be made in writing within 90 days of the date the product was shipped. Any claims not made within this time limit shall be deemed waived by the customer. Advancelab is not responsible for any additional costs of repair caused by poor packaging or in-shipment damage during return.

**Warranty Returns:** All warranty returns must be authorized in advance by Advancelab and approved by writing. Unless approved in advance for good reason, all returns must be in original condition, including all manuals, and must be packaged in original packaging materials. All returned goods are to be shipped to Advancelab, freight prepaid at customer's expense.



# WALK IN

AdvanceLab has built a rock steady reputation for solutions where most fail. Since our establishment in 2003, we have positioned ourselves as the premier facility solution provider for the scientific industry; designing, building and delivering quickly and without fuss. We have also taken our brand of expertise globally; now with exports of laboratory casework, fume hoods, laminar flow cabinets and clean booths heading to 30 countries, covering all continents.



# FOR YOUR FUTURE

## Indonesia

PT. Advancelab Saintifik

253 Kaki Bukit Ave 1  
Singapore 416619  
Tel: +65 6448 8255  
Fax: +65 6448 9833  
Email: [info@advancelab.com.sg](mailto:info@advancelab.com.sg)

## Malaysia

Advancelab Sdn Bhd

No. 3388, Jalan Pekeliling Tanjung  
27/2, Kawasan Perindustrian Indahpura,  
81000 Kulaijaya, Johor, Malaysia.  
Tel: +607 660 8877  
Fax: +607 660 8866  
Email: [info-my@advancelab-global.com](mailto:info-my@advancelab-global.com)

## Myanmar

Advancelab Scientific & Engineering Co., Ltd

No.(81/1)-1A, Myin Thar 7th Street,  
(14/1)Ward, South Okkalapa Township,  
Yangon, Myanmar  
Tel: +95 (9) 779753802  
Email: [info-mm@advancelab-global.com](mailto:info-mm@advancelab-global.com)

## Singapore

Advancelab (S) Pte Ltd

253 Kaki Bukit Ave 1  
Singapore 416619  
Tel: +65 6448 8255  
Fax: +65 6448 9833  
Email: [info@advancelab.com.sg](mailto:info@advancelab.com.sg)

## Thailand

Prima Scientific Co., Ltd

147 / 170 - 171 Baromrajchonnane  
Road, Arunamarin Bangkoknoi, Bangkok  
10700, Thailand.  
Tel: +66 2 884 9480  
Fax: +66 2 884 6441  
Email: [primasci@primasci.com](mailto:primasci@primasci.com)

Website:

[www.advancelabglobal.com](http://www.advancelabglobal.com)

Indonesia • Malaysia • Myanmar • Philippines • Singapore • Thailand

